



ELSEVIER

Contents lists available at ScienceDirect

## CYTOTHERAPY

journal homepage: [www.isct-cytotherapy.org](http://www.isct-cytotherapy.org)
 International Society  
**ISCT**  
 Cell & Gene Therapy®

## Review Article

## Manufacturing Processes

# A comprehensive review of quantum bioreactor cell manufacture: Research and clinical applications

 Charlotte H. Hulme<sup>1,2,\*</sup>, Claire Mennan<sup>1,2,\*</sup>, Helen S. McCarthy<sup>1,2</sup>, Rebecca Davies<sup>1,2</sup>,  
 Tian Lan<sup>1,2</sup>, Larissa Rix<sup>1,2</sup>, Jade Perry<sup>1,2</sup>, Karina Wright<sup>1,2,\*\*</sup>
<sup>1</sup> Centre for Regenerative Medicine Research, School of Pharmacy and Bioengineering, Keele University, Keele, Newcastle, United Kingdom

<sup>2</sup> Robert Jones and Agnes Hunt Orthopaedic Hospital, Gobowen, Oswestry, Shropshire, United Kingdom

## ARTICLE INFO

## Article History:

Received 19 December 2022

Accepted 10 April 2023

Available online xxx

## Keywords:

bioreactor

cell therapy

extracellular vesicles

Good Manufacturing Practice

Quantum

up-scale manufacture

## ABSTRACT

The Quantum cell expansion system manufactured by Terumo-BCT is perhaps the most widely reported Good Manufacturing Practice-compliant bioreactor used for the expansion of adherent cell populations, both for research purposes and clinical cell-based therapies/trials. Although the system was originally designed for adherent cell expansion, more recently suspension cultures and extracellular vesicle manufacturing protocols have been published using the Quantum system. Cell therapy research and regenerative medicine in general is a rapidly expanding field and as such it is likely that the use of this system will become even more widespread and perhaps mandatory, for both research and development and in the clinic. The purpose of this review is to describe, compare and discuss the diverse range of research and clinical applications currently using the Quantum system, which to our knowledge has not previously been reviewed. In addition, current and future challenges will also be discussed.

© 2023 International Society for Cell & Gene Therapy. Published by Elsevier Inc. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>)

## Introduction

The Quantum (Terumo-BCT, Lakewood, CO) is perhaps one of the most widely used bioreactor systems, to date, that is compliant for Good Manufacturing Practice (GMP). Here, we review the current status of the Quantum cell expansion system and give a flavor of its applications for research and clinical use. Further, we suggest strengths and weaknesses of the system from a user outlook and highlight future perspectives. Throughout, it should be noted that there may be other bioreactor options available for each of the applications that were outside the scope of this review.

## The Quantum Hollow-Fiber GMP Bioreactor

The Quantum cell expansion system is a functionally closed, automated, scalable, hollow-fiber bioreactor [1,2]. The bioreactor component comprises up to 11 500 hollow fibers, which create an internal surface area of 2.1 m<sup>2</sup>, equivalent to 120 T175 tissue culture flasks

[1,3]. Temperature/gas controls and critical processes such as feeding, removing waste products and harvesting cells are automated (Figure 1). The bioreactor has proven effective for the expansion of a variety of adherent cells, suspension cells, extracellular vesicles (EVs) and viral vectors (as in Table 1 [4–34]; discussed later).


### Why is Up-Scale Advanced Therapy Medicinal Products (ATMPs) Manufacture Needed?

Therapeutic doses of cell-based therapies can be very high depending on the application, for example, intravenous therapies often use >10<sup>9</sup> cells, requiring a scaled-up manufacture to produce adequate cell numbers and to ultimately reduce costs [4,35–37]. Many ATMP developers and manufacturers are moving away from autologous approaches to treatment and bespoke manufacture, in favor of allogeneic “off-the-shelf” approaches [38]. Up-scale manufacturing lends itself well to the need for large allogeneic banks of cells. Traditional biopharma facilities often use automated, closed technologies, which are designed to repeatedly manufacture identical products at high volume while balancing operation efficiency, cell doubling considerations, process stability and the risk of batch losses [36]. Flexibility is required as more ATMP products move toward gaining regulatory approval; however, cell therapies must be

\*\* Correspondence: Karina Wright, Centre for Regenerative Medicine Research, School of Pharmacy and Bioengineering, Keele University, Keele, Newcastle, ST5 5BG, United Kingdom.

E-mail address: [karina.wright1@nhs.net](mailto:karina.wright1@nhs.net) (K. Wright).

\* These authors contributed equally to this work.

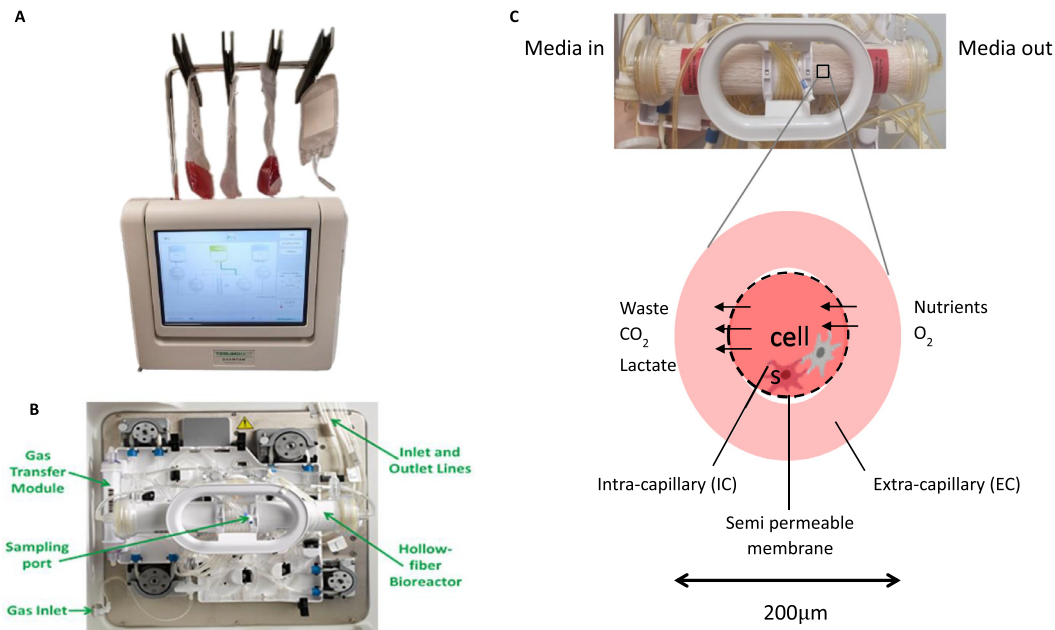

787

**Figure 1.** (A) The Quantum Cell Expansion System. (B) A single-use cassette, including the hollow-fiber bioreactor is inserted inside the Quantum Cell Expansion System, with the reactor connected to a gas supply and tubing inserted into the peristaltic pump system to move culture media through the system. (C) Schematic depicting the hollow fiber perfusion system inside the bioreactor, indicating the flow of nutrients, waste, CO<sub>2</sub>, O<sub>2</sub> and lactate across the semi-permeable membrane of the fibers.

**Table 1**  
Overview of Quantum-expanded cultures and substrates.

Species	Cell type	Number of studies	References
Human	BM-MSCs	14	[1,3,5,6,9–17,34]
	AD-MSCs	6	[18–23]
	UC-MSCs	2	[11,24]
	Neural MSCs	2	[4,25]
	Periosteum MSCs	1	[26]
	CD3+ T cells	2	[31,32]
	CD4+ T cells	1	[33]
	CD25+ T cells	1	[33]
	Dendritic cells	1	[29]
	Fibroblasts	1	[28]
	Skeletal muscle fibroblasts	1	[28]
	HEK 293T	1	[27]
	Embryonic stem cells	1	[2]
	IPSC endothelial cells	1	[8]
	IPSC	1	[7]
	Porcine	BM-MSCs	1
Murine	Embryonic fibroblasts	1	[2]
Substrates	Fibronectin	22	[1–6,8,10,13–17,19–22,24,26–28,30]
	Cryoprecipitate	3	[6,11,18]
	Vitronectin	2	[6,7]
	Laminin and vitronectin	1	[7]
Media	Alpha MEM	10	[1,6,10,12,14,15,18,20,24,30]
	DMEM	11	[2–5,11,13,16,17,21,22,26,27]
	Other custom XFSM	2	[7,28]
	PRIME-XV T-cell Expansion XFSM	2	[31,33]
	RPMI	1	[29]

AD-MSCs, adipose-derived mesenchymal stromal cells; BM-MSCs, bone marrow-derived mesenchymal stromal cells; DMEM, Dulbecco's modified Eagle's medium; HEK 293T, human embryonic kidney 293 SV40 T-antigen cells; IPSC, induced pluripotent stem cell; MEM, modified Eagle's medium; RPMI, Roswell Park Memorial Institute; UC-MSC, umbilical cord-mesenchymal stromal cells, XFSM, xeno- and serum free media.

produced using standardized, reproducible, GMP-compliant processes that can be tracked, documented and audited [39].

### Controlling and Monitoring Manufacture in the Quantum System

#### *Hypoxia and normoxia*

The physiologic oxygen concentration experienced by cells in the tissues of the body is markedly lower than the 21% oxygen tension commonly used for cell expansion *in vitro*. Previous research in our group and by others has demonstrated that the culture of cells in hypoxic conditions can be therapeutically beneficial in terms of retention of a desirable *in vivo* tissue-forming phenotype and a slowing of plastic adherence de-differentiation [40–44]. The Quantum bioreactor has a direct gas connection, which allows users to connect any blend of gas. Terumo-BCT states that an internal investigation demonstrated that a comparison of hypoxic and normoxic conditions resulted in no significant difference in cell yield, viability or phenotypic, morphological or differentiation potential, but the cell types grown were not clearly defined (personal communication, Terumo-BCT, October 2022). Hanley *et al.* [5] showed that the yield of bone marrow-derived mesenchymal stromal cells (BM-MSCs) was 75% greater in 5% oxygen cf. 21% oxygen in sister populations grown on tissue culture plastic. Further work needs to be undertaken to establish the use of maintaining a more “physiologic” oxygen tension during cell manufacture in the Quantum system.

#### *Coating the hollow fibers*

Cell adhesion is critical for the proliferation and differentiation of anchorage-dependent cells and relies on the surface chemistry of the tissue culture substrate. The hollow fibers of the Quantum are non-adherent and require coating with a substrate that is permissive to cell attachment and allows the cells to retain an *in vivo* phenotype. As summarized in Table 1, fibronectin, vitronectin and human donor pooled cryoprecipitate (CPPT) have most commonly been used as substrate coatings for the Quantum for the expansion of MSCs, whereas laminin or fibronectin have been used for the expansion of human-induced pluripotent stem cells [6–8]. Gelatin, human serum albumin, collagen I, poly-l-lysine and poly-d-lysine were also tested for MSCs but produced lower yields compared with the more commonly used substrates [6].

Fibronectin is a high molecular weight (440 kDa) glycoprotein found in the blood and as a solid component linking proteins in the extracellular matrix [45]. Plasma-derived fibronectin is widely used in cell culture, although recombinant and cellular fibronectin is also commercially available. Fibronectin benefits from being a fully defined protein that is suitable for multiple cell types. Vitronectin is a low molecular weight (75 kDa) serum protein that similarly interacts with cells via integrin-mediated adhesion to RGD motifs [46]. Currently, recombinant vitronectin and fibronectin products and other similar synthetic surface-coating reagents are commercially available that adhere to GMP manufacturing guidelines. However, when considering the use of fibronectin, the choice of commercial supplier requires careful consideration, with several vendors no longer being recommended for ATMP manufacture by the US Food and Drug Agency or requiring extensive, costly adventitious viral testing.

Cryoprecipitate is derived from human blood by freezing fresh plasma, then thawing and centrifuging to collect a precipitate. It is rich in clotting factors such as fibrinogen, factor VIII, factor XIII and von Willebrand factor, and it also contains fibronectin, platelet micro-particles and immunoglobulins. The fibrinogen concentration in CPPT varies from 10 g/L to 25 g/L, depending on the donor and the processing [47]. When CPPT is used as a Quantum substrate, it is usually prepared as a pooled product from four to six donors [47]. In summary, the efficiency of CPPT in coating the Quantum bioreactor differs

between batches, and the product is also likely to be more difficult to obtain than the other commercial bioreactor substrates described. Frank *et al.* [6] compared BM-MSCs that were Quantum expanded on fibronectin, recombinant vitronectin and CPPT. All substrates performed well for BM-MSC manufacture, with similar cell yield and no differences in growth kinetics, morphology, adherence to the International Society for Cell and Gene Therapy (ISCT) flow profile and capacity to trilineage differentiate. Therefore, the choice of coating reagent should be determined by the user dependent of their specific cell type, ideally tested on plastic before Quantum expansion and further influenced by whether the cell product needs to be manufactured in accordance with GMP.

#### *Media, monitoring cell growth and harvesting*

A plethora of culture media types have been used in the Quantum and have been established for different applications (summarized in Table 1). For example, BM-MSCs have been cultured in alpha-minimum essential media, Dulbecco's minimum essential media and GlutaMax (Gibco supplement, used as an alternative to L-glutamine), adipose-derived MSCs (AD-MSCs) in alpha-minimum essential media and GlutaMax, T cells in RPMI and a chemically defined xeno-free media (PRIME-XV T cell expansion XFSM media). Most types of media are compatible with the Quantum and are chosen to suit the cell type being cultured rather than the system itself. Certain cell types require additional supplements in addition to basal media to support growth in the Quantum and to induce or retain a desirable phenotype. These include fetal bovine serum (FBS), human platelet lysate (hPL), insulin—transferrin—selenium, glucose, vitamins (e.g., vitamin C [ascorbic acid] for chondrocytes) and amino acids (such as L-glutamine) [48]. The Quantum is designed to allow for the addition of cytokines, growth factors and supplements before and during the culture process enabling bespoke manufacturing processes.

Monitoring of cell growth kinetics in the system is commonly undertaken using lactate and/or glucose readings, sampled in the conditioned media acquired during a manufacturing growth phase. Sampling conditioned media via the sampling port for this purpose allows for analysis of the media currently circulating in the EC circulation loop. There are established protocols for using these metabolic measures (glucose and lactate) and that consider these alongside the perfusion rate through the bioreactor, in order to establish the cellular production rate per day [3]. This allows for prediction of local doubling time and, once worked up with matched tissue culture plastic (TCP) populations on the user's specific cells, can be used to predict the number of cells in the bioreactor system at any specific time. Plotting the calculated lactate production rate over time can allow the user to determine when the cells are in an expansion phase and, once this begins to plateau, can inform on an appropriate time to harvest. An additional/alternative method of assessing cell growth, indirectly, is to establish and monitor a “sister” population seeded onto traditional tissue culture plastic vessels (flasks/plates). Although not directly comparable with the conditions in the Quantum system (particularly in terms of substrate and media perfusion), such parallel cultures allow the operator to evaluate how well a particular culture is growing at the low densities most commonly used in the bioreactor and can inform critical decisions such as need to alter perfusion rates and time to harvest. Harvesting cultures from the Quantum is usually undertaken via standard trypsinization for adherent cell populations with the use of animal derived trypsin or recombinant TrypLE (Thermo Fisher Scientific, Waltham, MA, USA). A range of concentrations have been used previously; however, 0.25% solutions tend to be preferred for adherent bioreactor cell cultures cf. the 0.05% concentration more commonly used in traditional tissue culture plastics.

## Cell and EV Manufacture

### *Applications of Quantum cell manufacture for research*

The Quantum system has been widely used in the up-scale expansion of many different cell types for research purposes. Some have been aimed at yielding maximal numbers of cells from a single passage, whereas others have looked to the Quantum as a method of automating processes in developing protocols that will be transferable for GMP cell manufacture. Most of the research studies reviewed in this manuscript have focused on the derivation and expansion of BM-MSCs or from subcutaneous fat/adipose (AD-MSCs), isolated and/or expanded on fibronectin (Table 1). For BM-MSC isolation, BM aspirate can be loaded directly into the Quantum and large numbers of BM-MSCs (in the region to 10–35 million cells) yielded from the first passage [1,5,9–17]. Once isolated, these BM-MSCs have been reseeded into the Quantum and further culture expanded, with cell yields typically being greater than 100 million MSCs [1,5,9,11,17].

Several studies have demonstrated that Quantum-derived BM-MSCs adhere to the ISCT guidelines in terms of cell-surface marker profile and trilineage differentiation potential [1,5,9,11,14–16,49,50]. Functionality of Quantum-manufactured BM-MSCs has been demonstrated using *in vitro* activated T-cell suppression assays [5,12,15,16]. Further, the efficacy of Quantum-derived BM-MSCs has been demonstrated *in vivo* in rodent models for the treatment of ischemic stroke and cartilage joint surface injury [5,51].

Important considerations for GMP manufacturing process development and safety are whether Quantum expansion is successful using cryopreserved banks and determining whether expansion in the system induces any chromosomal aberrations. Rojewski *et al.* [1] demonstrated that BM-MSCs which had previously been cryopreserved could be expanded in the Quantum, however their doubling time was 1.4-fold higher in comparison to seeding freshly isolated cells. In addition, Jones *et al.* [10] and Mennan *et al.* [11] demonstrated no detectable changes in Quantum-manufactured cultures *cf.* BM-MSCs expanded in parallel on standard TCP in terms of telomere length and spectral karyotype, micronucleus formation and tumorigenicity.

Process optimization in the Quantum is continually being sought in terms of limiting the frequency of “open” procedures, improved automation and reducing costs. Hanley *et al.* [5] calculated that for a theoretical clinical trial using BM-MSCs to treat 60 patients, the use of the Quantum required 0.02% of the number of “open” procedures compared with TCP, thus significantly reducing the likelihood of product contamination. This study also indicated that, at the time of publication, the cost of manufacturing 200 million BM-MSCs was comparable in either the Quantum or TCP; however, the labor costs were much lower, primarily due to the reduction in time to harvesting the cells from the bioreactor (approx. 45 min), in comparison with TCP (approx. 6 hours) [5]. Further, Russell *et al.* [15] indicated that the manufacture of 100 doses of 100 million BM-MSCs would be manufactured a month earlier using Quantum *cf.* TCP and that the reagents and consumables costs would be nearly halved by using Quantum, equating to a saving of \$976.41 per dose.

In addition to BM-MSCs, MSCs derived from alternative tissues have been widely expanded in the Quantum system (Table 1). Following BM, the most frequently used source of MSCs that have been expanded in the bioreactor are AD-MSCs [18–23]. These studies have included loading raw liposuction aspirate into the Quantum, isolating and expanding the AD-MSCs in the system [18–20] as well as re-seeding Quantum-isolated AD-MSCs [18,20] and seeding AD-MSCs prior isolated and expanded on TCP [21,22]. MSCs in the order of hundreds of millions are typically harvested following Quantum isolation, with reported cell viabilities of >95% [18,21]. Adherence to the ISCT MSC profile was maintained following AD-MSC Quantum expansion [18–22]. Other cell types successfully expanded in the

Quantum system include MSCs isolated from umbilical cord digests [11], Wharton's jelly [24], neural tissues [25] and periosteum [26], human iPSCs [8], HEK293T cells [27], fibroblasts and skeletal myoblasts [28] monocyte-derived dendritic cells [29], as well as animal-derived cells such as mouse embryonic fibroblasts [2] and porcine MSCs [30].

Besides MSCs, there are also several reports of T-cell expansion protocols in the Quantum [31–33]. We realize that the expansion of T-cell populations is far more established in other bioreactor systems, for example, the Xuri Cell Expansion system (Cytiva, Marlborough, MA, USA) [52], the CliniMacs Prodigy (Miltenyi Biotec, Bergisch Gladbach, Germany) [53,54] and the Cocoon cell therapy manufacturing platform (Lonza, Basel, Switzerland) [55] and the use of the Quantum for this purpose is still in early process development. However, we thought it worthy to include a few examples of this alternate use of the system in our review. In brief, Nankervis *et al.* [32] optimized the Quantum expansion process for CD3+ T-cells, demonstrating that perfusion feeding allowed efficient gas exchange and continued feeding of cells which in turn increased yields. Coeshott *et al.* [31] further demonstrated that CD3+ T cells could be expanded from the Quantum when seeded at low ( $30 \times 10^6$  cells) or high ( $85 \times 10^6$  cells) densities, yielding  $1.48 \times 10^{10}$  and  $2.51 \times 10^{10}$  cells, respectively. The 8–9 days required to culture expand these large yields has the potential to reduce the production time for CAR-T cell therapies and indicates that lower seeding densities could be used for individuals with low numbers of circulating CD3+ T cells [31]. It has also recently been shown that the Quantum cells can more efficiently support the up-scale expansion of CD4+ and CD25+ regulatory T cells in comparison with standard TCP [33].

### *Applications of Quantum for EV manufacture*

So far, we have reviewed the multiple and diverse cell expansion protocols using the Quantum system, but an emerging additional purpose for the Quantum is to enhance the production of cell-derived products such as EVs. EVs are small cell-derived particles (<1  $\mu$ m) that contain mRNA, miRNA and protein cargo that is encapsulated by a lipid membrane [56]. Protocols have been developed to enrich EVs (particularly those <200 nm in size) for use as an alternative to cell-based therapies [57]. However, there are still hurdles to overcome in the manufacture of clinically applicable EVs.

The Quantum platform represents an opportunity to overcome current limitations of the field, mainly poor vesicle yields, of which the resulting EVs are highly heterogeneous, labor-intensive and vary between batches. Already, some studies have begun to establish this reality, being able to scale-up the production of EVs, typically of MSC origin. These studies differ by many factors that have yet to be standardized, including different harvest medias (either EV-depleted FBS/hPL or serum free to minimize contaminating EVs) and isolation protocols (Table 2) [12,24,58–61]. Although this is a common caveat in EV research, the Quantum requires further research to devise an optimal protocol to harvest EVs. Such factors for exploration should include the optimal time-point to harvest EVs, informed by glucose and lactate readings, and the chosen strategy for collection, which can be done continuously, in batches, or at the end of the run. As more studies emerge, full disclosure of the used materials, protocol and experimental outputs (metabolic measures, predicted cell number and cell number at harvest) in consideration to EV outputs will be integral.

What remains to be established is whether Quantum-derived EVs retain their therapeutic attributes and how this changes the final product, since culture conditions can greatly affect EV characteristics [62]. Already, Quantum MSC EVs have been shown to significantly improve neurological recovery in swine or monkey models of traumatic brain injury [58,60,61]. However, such studies based this conclusion in comparison with vehicle controls, and no direct

**Table 2**  
Extracellular vesicle manufacture in the Quantum.

Cells loaded for expansion	Expansion media	EV harvest media	EV harvest frequency	EV isolation protocol	EV characterization	EV characteristics	Reference
$3 \times 10^6$ human BM-MSC	N/D		Start on day 3, followed by every 48 h for 6 d and every 24 h for the final 4 d	250g for 5 min → 3000g for 30 min → 0.22- $\mu$ m filtration → 100 000g for 120 min	qNano (IZON)	N/D	[58]
$20 \times 10^6$ human BM-MSC	$\alpha$ MEM + 1% L-glutamine + 5% HPL	$\alpha$ MEM + 1% L-glutamine	Every 48 h for 12 d once peak expansion was achieved	1000g for 15 min → 0.2- $\mu$ m filtration → 100 000g for 180 min (Ti45 rotor)	NanoSight; Protein assay; Flow cytometry	$9.8\text{--}44.9 \times 10^{11}$ per harvest (0.9–6.2 mg protein); $9.8\text{--}15.6 \times 10^{13}$ total particle yield, CD9+, CD63+, CD81+ and CD47+	[59]
$10 \times 10^6$ monkey BM-MSC	$\alpha$ MEM + 10% EV-depleted FBS		Every 48 h for 4 d, then every 24 h for 2 d	250g for 5 min → 3000g for 30 min → 0.22- $\mu$ m filtration → 100 000g for 120 min	qNano (IZON)	111-nm mode size	[60]
$10 \times 10^6$ monkey BM-MSC	$\alpha$ MEM + 10% EV-depleted FBS		Start on day 3, followed by every 48 h for 4 days and every 24 h for the final 2 d	250g for 5 min → 3000g for 30 min → 0.22- $\mu$ m filtration → 100 000g for 120 min	qNano (IZON)	N/D	[61]
$30 \times 10^6$ human WJ-MSC	$\alpha$ MEM + [ 20% FBS OR 5% hPL]	$\alpha$ MEM	Every 24 h for 3 d once peak expansion was achieved	2000g for 30 min → 20 000g for 30 min → 100 000g for 70 min (Ti70 rotor)	NanoSight; Flow cytometry	$1.13 \pm 0.5 \times 10^{13}$ , [FBS] $1.64 \pm 0.2 \times 10^{13}$ [hPL] total particle yield, $113 \pm 6.2$ nm [FBS], $125 \pm 4.4$ nm [hPL], mode size, CD63+ and CD81+	[24]
$10 \times 10^6$ monkey BM-MSC	$\alpha$ MEM + 10% EV-depleted FBS		Start on day 3, followed by every 48 h for 4 d and every 24 h for the final 2 d	250g for 5 min → 3000g for 30 min → 0.22- $\mu$ m filtration → 100 000g for 120 min (JS-24.38 swing bucket rotor)	qNano (IZON); TEM; Western blot	30- to 150-nm size range, peaking at 100 nm, CD63+ and Alix+	[12]

$\alpha$ MEM, alpha minimum essential media; BM-MSC, bone marrow–derived mesenchymal stromal cells; EV, extracellular vesicle; FBS, fetal bovine serum; hPL, human platelet lysate; WJ-MSC, Wharton's jelly–mesenchymal stromal cells.

comparison has been made to tissue culture plastic, for example. These studies also lack adequate EV characterization in accordance with the International Society of Extracellular Vesicles minimal reporting criteria [63,64]. Such criteria are necessary since, to date, no existing protocol can isolate pure EVs, and it is important to establish what co-isolates may contribute to findings. Zhang *et al.* [12] have better accommodated this, associating CD63+ and Alix+ Quantum-derived particles of 30–150 nm with improved cognitive function in the same traumatic brain injury model, but further characterization is necessary to fully define a Quantum EV preparation. Here, following International Society of Extracellular Vesicles recommendations should dictate the minimum information on EVs, but future research should focus on a wider-view by using 'omic techniques. This will be crucial in moving towards the use of the Quantum system to manufacture EVs for clinical application.

In this regard, several aspects of EV manufacture will need to be adjusted for GMP compliance, and the Quantum system represents an opportunity to move the field further toward a clinically translational product. Bellio *et al.* [24] have begun some of the work necessary, exploring the difference between Wharton's Jelly–MSCs expanded in 20% FBS or 5% hPL before EV harvest in serum-free media, demonstrating no significant differences in particle yield ( $1.13 \times 10^{13}$  particles/mL, FBS versus  $1.64 \times 10^{13}$  particles/mL, hPL) or mode size (113 nm FBS versus 125 nm hPL), and both preparations were shown to express EV-specific markers, CD63 and CD81, with few differences being found in their miRNA cargo. Furthermore, both EV enrichments improved cardiac function post-myocardial infarction. Mendt *et al.* [59] have taken this work a step further, using the Quantum to achieve GMP-grade EVs, by transfer of conditioned media from the bioreactor in a closed system to conduct isolation,

in accompaniment to sterility testing. This generated a total of  $9.8\text{--}15.6 \times 10^{12}$  particles per run, with  $9.8\text{--}44.9 \times 10^{11}$  particles being collected every 48 hours depending on the day of harvest. Here, EVs were generated for the purpose of making "iExosomes," engineered EVs with the ability to target oncogenic Kras, used to treat patients with pancreatic cancer [59]. Here, both the reporting of Quantum outputs (matching particle counts to lactate readings [24] and characterization of individual batches [59]) and EV characterization is much improved. However, insights into protocol development, the effect on EVs due to their Quantum origin and consideration to potential non-EV particulate is still necessary to address. As such, whilst GMP grade EVs are achievable, these avenues would need to be explored to make the Quantum an established source of clinically applicable EVs.

#### Clinical applications of Quantum-derived cells

We have independently searched for publications and trials in which the Quantum system is used clinically and have also requested all available information from Terumo-BCT. Clinical trial information available on the NCT database was not always detailed enough to divulge the cell manufacture methodology, and therefore this summary is unlikely to be exhaustive. Here we present a flavor of the range of clinical applications currently undertaken in the system (Table 3) [4,9,23,34,65–74], but it is beyond the scope of this review to fully analyze the findings from each study. Clinical trials using the Quantum explore a range of cell types to evaluate their efficacy for varying indications. The use of a bioreactor such as the Quantum can enable more centers to treat patients with cell-based therapies, where an on-site GMP compliant laboratory is not available. One such example is a case report of a 24-year-old man with grade III graft-versus-host disease (GVHD) following an allogeneic stem cell

**Table 3**  
Quantum-manufactured cells in the clinic.

Cell type	Indication	Cell yield/dose	Period in Quantum	Trial ID	Reference
BM-MSCs	Graft-versus-host disease (GVHD)	Two doses of $1.5 \times 10^6$ cells/kg	23 d (2 passages)	N/A (hospital exemption)	[9]
	Immune modulation in patients receiving full-thickness corneal transplants	$2 \times 8.0 \times 10^7$ or $2 \times 1.6 \times 10^8$	6–8 d	Not registered	Ongoing trial, no results published to date (VISICORT)
	Type 2 diabetes and DKD	$8 \times 10^7$ , $1.6 \times 10^8$ or $2.4 \times 10^8$ fixed doses of ORBCELL-M	N/D	NCT02585622	Ongoing trial, no results published to date (NEPHSTROM)
	Idiopathic Parkinson disease	$1 \times 10^6$ cells/kg, $3 \times 10^6$ cells/kg, $6 \times 10^6$ cells/kg or $1 \times 10^7$ cells/kg	3 passages	NCT02611167	[34]
GM-neural MSCs with cytosine deaminase (CD-NSCs)	Phase 1 trial for recurrent high-grade gliomas	$1 \times 10^7$ or $5 \times 10^7$	7–10 days +9 days	NCT01172964	[4,65]
	Phase 1 trial for recurrent high-grade gliomas	$5 \times 10^7$ to $1.5 \times 10^8$	7–10 days +9 days	NCT02015819	[66]
GM-neural MSCs with carboxylesterase (CE-NSCs)	Phase 1 trial for recurrent high-grade gliomas	N/D	N/D	NCT02192359	Ongoing trial, no results published to date
AD-MSCs	Allogeneic adipose-derived stromal cells for patients with ischemic heart disease and ischemic heart failure (phase 1)	$1 \times 10^8$ cells	2 passages	NCT02387723	[67]
	Allogeneic adipose-derived stromal cells for patients with ischemic heart disease and ischemic heart failure (phase 2)	$1 \times 10^8$ cells	2 passages	NCT03092284	[68]
	Stem Cell Therapy in Ischemic Non-treatable Cardiac Disease (SCIENCE)	$1 \times 10^8$ cells	2 passages	NCT02673164	[69]
	Primary graft dysfunction after lung transplantation	$1 \times 10^8$ cells	2 passages	NCT04714801	Ongoing trial, no results published to date
	Aqueous deficient dry eye disease	$22 \times 10^6$ MSCs/mL	N/D	NCT03878628	[23]
	MultiStem (BM-derived Multipotent Adult Progenitor Cells, or MAPC)	MultiStem in Ischemic Stroke (MASTERS)	$4 \times 10^8$ or $1.2 \times 10^9$ cells	16 d (2 passages)	NCT01436487
MultiStem Administration for Stroke Treatment and Enhanced Recovery Study (MASTERS-2)		$1.2 \times 10^9$ cells	16 d (2 passages)	NCT03545607	Ongoing trial, no results published to date
Treatment Evaluation of Acute Stroke for Using in Regenerative Cell Elements (TREASURE)		$1.2 \times 10^9$ cells	16 d (2 passages)	NCT02961504	Rationale in [72] Ongoing trial, no results published to date
Efficacy and Safety Study of HLCM051 (MultiStem) for Pneumonic Acute Respiratory Distress Syndrome (ONE-BRIDGE)		$9 \times 10^8$ cells	16 d (2 passages)	NCT03807804	[73]
Phase 1/2 Study to Assess MultiStem Therapy in Acute Respiratory Distress Syndrome (MUST-ARDS)		$9 \times 10^8$ cells	16 d (2 passages)	NCT02611609	[74]
Fibroblasts and keratinocytes	Autologous dermo-epidermal skin substitute for the treatment of burns	N/D	N/D	NCT03227146, NCT03229564,	Ongoing trial, no results published to date
	Autologous dermo-epidermal skin substitute for the treatment of full-thickness skin defects	N/D	N/D	NCT03394612	Ongoing trial, no results published to date

BM-MSC, bone marrow–derived mesenchymal stromal cells; DKD, diabetic kidney disease; N/A, not available; N/D, not determined.

transplant [9]. Because of the severity of the condition and the speed at which treatment was required, BM-MSCs were grown in the Quantum bioreactor. A 25-mL sample of BM aspirated from the original donor yielded  $6 \times 10^6$  cells following 23 days in culture (2 passages), and two infusions of  $1.5 \times 10^6$  cells/kg was found to reduce the severity of GVHD to grade I in 3 months.

VISICORT is evaluating the potential for Quantum allogeneic BM-MSCs to modify the immune system's response for patients receiving a full-thickness corneal transplant. In this phase 1b trial, healthy BM-MSCs are expanded in the Quantum for 6–8 days, cryogenically preserved in individual dose infusion bags. This study reports on eight patients, in whom the safety and tolerability of two intravenous

infusions of cryopreserved BM-MSCs at 7 days and 1 day before the corneal transplant procedure are tested. Group 1 ( $n = 4$ ) will receive two infusions of  $8 \times 10^7$  cells and group 2 ( $n = 4$ ) will receive  $1.6 \times 10^8$  cells. Subject to satisfactory results, a larger phase 2 trial is planned.

The Novel Stromal Cell Therapy for Diabetic Kidney Disease (NEPHSTROM) study is a double-blind, placebo-controlled phase 1/2 trial investigating the safety, feasibility, tolerability and preliminary efficacy of an intravenous infusion of ORBCELL-M (Orbsen Therapeutics, Galway, Ireland) Quantum manufactured BM-MSCs (<http://clinicaltrials.gov: NCT02585622>). Up to 48 patients with type 2 diabetes and progressive kidney disease will be exposed to either placebo,  $8 \times 10^7$ ,  $1.6 \times 10^8$  or  $2.4 \times 10^8$  BM-MSCs ( $n = 12$  each group), and markers of diseases assessed for up to 18 months' post-infusion. Another phase 1 study assessed the safety and tolerability of intravenous allogeneic BM-MSCs in patients with idiopathic Parkinson disease (<http://clinicaltrials.gov: NCT02611167>) [34]. Here, 20 patients were sequentially enrolled and received a single infusion of either  $1 \times 10^6$  cells/kg,  $3 \times 10^6$  cells/kg,  $6 \times 10^6$  cells/kg or  $1 \times 10^7$  cells/kg. The study concluded infusions at these concentrations were safe, well tolerated and not adversely immunogenic [34].

Allogeneic neural stem cells (NSCs) have been expanded in the Quantum to produce clinical-grade genetically modified NSCs stably expressing the prodrug-activating enzyme cytosine deaminase (CD-NSCs), and subsequent adeno-viral transduction to express a modified human carboxylesterase (CE-NSCs) [4]. A single intracerebral administration of either  $1 \times 10^7$  or  $5 \times 10^7$  CD-NSCs in 15 patients with recurrent high-grade gliomas in a phase 1 trial has successfully demonstrated safety, migratory ability of the cells to tumor sites and successful conversion of a subsequent orally administered prodrug 5-fluorocytosine into an active chemotherapy agent (<http://clinicaltrials.gov: NCT01172964>) [65]. Consequently, a further phase 1 trial investigated the feasibility of multiple doses of CD-NSCs; 15 patients received a median of 4 (range 2–10) intracerebral CD-NSC doses, ranging from  $5 \times 10^7$  to  $1.5 \times 10^8$  CD-NSCs in combination with 5-fluorocytosine  $\pm$  leucovorin, with no clinical signs of immunogenicity (<http://clinicaltrials.gov: NCT02015819>) [66]. The authors recommended a dose of  $1.5 \times 10^8$  CD-NSCs to be used in phase 2 trials. A further, ongoing phase 1 trial due for completion mid-2023, is investigating the safety and therapeutic dose of intracranially administered CE-NSCs (<http://clinicaltrials.gov: NCT02192359>).

The Cardiology Stem Cell Centre in Rigshospitalet, University Hospital Copenhagen, Denmark, holds a Danish manufacturing authorization (no. 23909) and a tissue establishment authorization (no. 32298) and produce GMP-grade allogeneic Quantum expanded ADMSCs (certificate of GMP compliance no. DK IMP 92217). These ADMSCs have recently been used in both phase 1 and 2 trials in patients with ischemic heart disease and ischemic heart failure, (<http://clinicaltrials.gov: NCT02387723>, [NCT03092284](http://clinicaltrials.gov: NCT03092284) and [NCT02673164](http://clinicaltrials.gov: NCT02673164)). In each case,  $1 \times 10^8$  cryopreserved allogeneic AD-MSCs were intramyocardially injected into the peri-ischemic area, in a series of small injections (12–20) with the aim of regenerating the damaged myocardium; initial safety and feasibility studies demonstrated no complications or serious adverse events [67]. Results from two randomized double-blind placebo-controlled phase 2 trials, one single-center (81 patients) and one Horizon 2020-funded multi-center European trial (138 patients) for the same clinical indication, are yet to be published (for trial design, see [68,69]). Manufactured AD-MSCs from the Cardiology Stem Cell Centre have been used at the same dose ( $1 \times 10^8$  cells) in a phase 1/2 double-blind placebo trial in 30 patients undergoing lung transplantation, with the aim of reducing post-transplantation immunological reactions and ischemic reperfusion injury (<http://clinicaltrials.gov: NCT04714801>) and at  $2.2 \times 10^7$  ASCs/mL in a phase 1 trial in 7 patients with aqueous deficient dry eye disease (<http://clinicaltrials.gov: NCT03878628>) [23].

MultiStem (Athersys, Inc., Cleveland, OH, USA) is a proprietary commercial pharmaceutical-grade composition of Multipotent Adult Progenitor Cells (MAPCs) specially formulated for clinical use by the clonal expansion of BM-MSCs in the Quantum. A successful randomized, double-blind, placebo-controlled phase 2 trial for patients suffering from ischemic stroke (MultiStem Administration for Stroke Treatment and Enhanced Recovery Study [MASTERS]; <http://clinicaltrials.gov: NCT01436487> [70,71]) demonstrated a safe intravenous administration of up to  $1.2 \times 10^9$  cells, 24–48 hours post-symptom onset. This has led to two further phase 3 trials (<http://clinicaltrials.gov: NCT03545607> and [NCT02961504](http://clinicaltrials.gov: NCT02961504)) to fully evaluate the efficacy of MAPC in the treatment of ischemic stroke [72]. MultiStem has also been used in a series of clinical trials for acute respiratory distress syndrome. Phase 1 and 2 trials across the United States and Japan for pneumonia- and coronavirus disease 2019-induced acute respiratory distress syndrome (<http://clinicaltrials.gov: NCT03807804> [73] and [NCT02611609](http://clinicaltrials.gov: NCT02611609) [74]) have demonstrated a tolerance for the intravenous administration of up to  $9 \times 10^8$  MAPCs and have led to a phase 2/3 trial with 400 patients, due for completion in December 2023 (<http://clinicaltrials.gov: NCT04367077>). denovoSkin is an autologous bio-engineered skin graft developed for use in severe burns and skin reconstructive surgery, whereby autologous fibroblasts and keratinocytes are expanded in the Quantum prior to being seeding onto a hydrogel scaffold (CUTISS AG, Schlieren, Switzerland) [75]. Phase 2 trials, which were due for completion by December 2022, are currently underway in adults and children for severe skin burns (<http://clinicaltrials.gov: NCT03227146> and [NCT03229564](http://clinicaltrials.gov: NCT03229564)) and full-thickness skin defects (<http://clinicaltrials.gov: NCT03394612>).

It is clear that the use of the Quantum bioreactor is gaining traction in providing adequate cell numbers for application in clinical trials. The safety and efficacy of the cells produced by the system has been demonstrated many times, indicating the system is suitable for producing large numbers of cells while maintaining phenotypic characteristics.

## Discussion

The up-scale expansion of MSCs using traditional TCP techniques remains an ongoing challenge. Large-scale cell manufacture requires a level of automation while retaining high cell quantity and quality. The surface area provided for cell attachment by the hollow fiber design of the Quantum is estimated to support a theoretical cell yield of  $1000 \times 10^6$  MSCs; however, cell yields achieved in the literature vary significantly from  $12 \times 10^6$  to  $605 \times 10^6$  MSCs [17,18]. Harvests for non-adherent cells in the Quantum are even greater, with reports of up to  $25 \times 10^6$  T cells [31]. Side-by-side comparisons of growth rates in sister populations grown on TCP have shown a faster growth rate in the Quantum [1,5,11,17,20,26]. Further, analyses of cells grown in TCP and the Quantum have demonstrated comparable characteristics in terms of phenotype and functionality [1,5,11,17,20,26]. The higher growth rate observed in the Quantum, however, could be influenced by the continuous perfusion of nutrients in the bioreactor cf. intermittent media changes in TCP.

Most studies using the Quantum have the end goal of developing a product that will be compliant with GMP manufacture. These studies have included manufacture for allogeneic and autologous applications, both of which have utility, such as when rapid, up-scale expansion of autologous cells is required, e.g., for GVHD and for manufacturing large numbers of cell doses for “off-the-shelf” allogeneic cell applications. Some current GMP manufacture includes the use of xeno- and blood products, for example, FBS. However, to future-proof these protocols, additional scrutiny and more rigorous testing may be needed and ideally protocols developed which move away from the use of xeno- and blood products. In order to coat the Quantum fibers, there needs to be a move away from human

**Table 4**  
Advantages and disadvantages of the Quantum bioreactor system.

	Advantages	Disadvantages
Cost	Reduced labor costs compared with traditional TCP expansion	High start-up costs and ancillary equipment/consumables
Single-use bioreactor	Time efficient No risk of cross-contamination Process flexibility	Expensive Plastic waste/disposal Not fully customizable
Hollow fiber system	Large surface area with continuous perfusion Single-user can operate multiple bioreactors	Cannot monitor cell morphology or visualize cell growth No online control to allow assessment offsite Requires relatively high starting cell yield Not fully optimized for Extracellular vesicles or suspension cultures Requires large media volumes held at room temperature
Regulatory	Reduced chance of contamination due to the closed system GMP-compliant; CE marked and FDA, MHRA compliant Reduces need for costly clean rooms	Difficult to detect a contamination early

CE, Conformance Européenne; FDA, Food and Drug Administration; GMP, Good Manufacturing Practice; MHRA, Medicines and Healthcare products Regulatory Agency; TCP, tissue culture plastic.

cryoprecipitate. Therefore, further studies similar to those performed by Frank *et al.* [6] are required to compare the efficiency of GMP substrates. Moreover, identification of serum-free media that can support cells within the hollow fiber system is required. The use of commercially available serum-free media will have significant cost implications for the manufacture of cells within the Quantum system, due to the large media volumes required. Therefore, in-house development of serum free media may be beneficial.

When optimizing protocols within the Quantum, the cost of the systems consumables and reagents can pose a significant challenge. There may be a reluctance to test multiple different optimization steps and technical replicates, due to the prohibitive cost of the bioreactors, their ancillary products (e.g., media bags) and the affiliated reagent costs of the large volumes of media, sera and coating reagents required. There are several ancillary products required for manufacture in the system, which include specialized bags for holding media/phosphate-buffered saline, waste conditioned media and for loading the cells/reagents. These ancillary products must be purchased through Terumo-BCT, thus maintaining non-competitive pricing and potentially limiting how “future-proofed” novel products may be should the company cease selling these products. At each step where media, reagents or phosphate-buffered saline must be added to the specialized ancillary bags, an “open” procedure must take place in which these reagents are “drawn into” the bags using a specialized peristaltic pump; potentially opening the system to contamination. These loaded bags are then “sterile welded” onto the disposable bioreactor cassettes. However, as the popularity of the system gains momentum for clinical applications, other biological companies (e.g., Thermo Fisher Scientific pre-culture filled cell media bags) are starting to develop products that could be directly loaded onto the Quantum in the future, including specialized GMP-compliant media pre-loaded into bags. Moreover, the starting minimum cell seeding density recommended for Quantum expansion can also represent a limiting factor, particularly when working with cells derived from primary tissues or rare specimens. Although hybrid processes can be used, in which the cells are isolated and culture expanded on tissue culture plastic before seeding within the bioreactor, it is not always desirable for cells to have undergone large numbers of preceding population doublings. A hybrid approach also increases the number of “open” procedures and moves away from optimal GMP manufacturing protocols. To address these issues, Terumo-BCT have now developed a new smaller bioreactor, called the Quantum Flex. The Quantum Flex system offers the opportunity to use fewer cells (<30 million T cells or <5 million MSCs), allowing for easier process development before switching to the larger sized bioreactors required for up-scale manufacturing. The utility of this system will become apparent as more studies using this product are published.

The Quantum bioreactor is popular and widely used across academia, in health care settings and industry. It is proposed that one of

the major advantages of the Quantum system is that it can significantly reduce costs associated with the required numbers of expert GMP manufacturing staff who would be required to maintain the large numbers of tissue plastic cultures comparable with the surface area of the Quantum bioreactor. In our experience as users, it is possible that a single user could manage the oversight of multiple Quantum systems simultaneously. However, a period of expert training is required to fully develop a working knowledge of the technology and how to troubleshoot both hardware and software errors in order to maintain optimal performance of the instrument. Additionally, we are aware that the use of multiple bioreactors reflects a “scale out,” rather than a “scale up” of manufacturing processes.

Although there are many advantages of this system, as discussed within this review, improvements could be made to further benefit the user and to ensure a more consistent manufactured product, as illustrated in Table 4. Currently, the Quantum does not have any monitoring system or alarm for the gas perfusion system. Hence, the user must ensure that the gas supply can be accessed via the input regulator and by monitoring the cells and media. As phenol red is not typically desirable for GMP manufacture, this may pose a difficulty (personal communication with ATMP manufacturing facilities, September 2022). Moreover, glucose and lactate measurements to determine cellular metabolic activity of the cells expanded in the Quantum are undertaken via daily conditioned media sampling from the outlet port and testing outside of the system using ancillary equipment. It would be beneficial for the Quantum to have an internalized means of testing and recording these parameters, such that the user could be automatically alerted to change the media flow rate in response to metabolic activity changes, or even better, that the system could automatically “self-monitor” the cells growth to ensure the cells are maintained in an active growth phase.

### Conclusions

This review of existing literature has highlighted the wide-scale use of the Quantum bioreactor in research and translational settings. We have highlighted a need to develop standardized protocols for each cell (or cell-derived) Quantum product, allowing for direct comparisons across studies/trials, and the need to ensure consistent, well-characterized endpoints, essential for the clinical application of Quantum-derived products.

### Author Contributions

Conception and design of the study: CH, CM, HSM, RD, TL, LR, JP and KW. Acquisition of data: CH, CM, HSM, RD, TL, LR, JP and KW. Analysis and interpretation of data: CH, CM, HSM, RD, TL, LR, JP and KW. Drafting or revising the manuscript: CH, CM, HSM, RD, TL, LR, JP and KW. All authors have approved the final article.



## Declaration of Competing Interest

The authors have no commercial, proprietary or financial interest in the products or companies described in this article.

## Funding

The authors thank the funders for their support of this work, including Versus Arthritis (grants 18480 and 21156), the Medical Research Council (grants MR/N02706X/1 and MR/S015167/1 and Confidence in Concept Award [MC\_PC\_16049]) and the Orthopaedic Institute, Oswestry (grant RPG 188) for supporting the salaries of the authors.

## References

- [1] Rojewski MT, Fekete N, Baila S, Nguyen K, Fürst D, Antwiler D, et al. GMP-compliant isolation and expansion of bone marrow-derived MSCs in the closed, automated device quantum cell expansion system. *Cell Transplant* 2013;22:1981–2000.
- [2] Roberts I, Baila S, Rice RB, Janssens ME, Nguyen K, Moens N, et al. Scale-up of human embryonic stem cell culture using a hollow fibre bioreactor. *Biotechnol Lett* 2012;34:2307–15.
- [3] Martin-Manso G, Hanley PJ. Using the quantum cell expansion system for the automated expansion of clinical-grade bone marrow-derived human mesenchymal stromal cells. *Methods Mol Biol* 2015;1283:53–63.
- [4] Tirughana R, Metz MZ, Li Z, Hall C, Hsu D, Beltzer J, et al. GMP production and scale-up of adherent neural stem cells with a quantum cell expansion system. *Mol Ther Methods Clin Dev* 2018;10:48–56.
- [5] Hanley PJ, Mei Z, Durett AG, da Graca Cabreira-Harrison M, Klis M, Li W, et al. Efficient manufacturing of therapeutic mesenchymal stromal cells with the use of the Quantum Cell Expansion System. *Cytotherapy* 2014;16:1048–58.
- [6] Frank ND, Jones ME, Vang B, Coeshott C. Evaluation of reagents used to coat the hollow-fiber bioreactor membrane of the Quantum<sup>®</sup> Cell Expansion System for the culture of human mesenchymal stem cells. *Mater Sci Eng C Mater Biol Appl* 2019;96:77–85.
- [7] Paccola Mesquita FC, Hochman-Mendez C, Morrissey J, Sampaio LC, Taylor DA. Laminin as a potent substrate for large-scale expansion of human induced pluripotent stem cells in a closed cell expansion system. *Stem Cells Int* 2019;2019:9704945.
- [8] Hochman-Mendez C, Mesquita FCP, Morrissey J, da Costa EC, Hulsman J, Tang-Quan K, et al. Restoring anatomical complexity of a left ventricle wall as a step toward bioengineering a human heart with human induced pluripotent stem cell-derived cardiac cells. *Acta Biomater* 2022;141:48–58.
- [9] Kozanoglu I, Maytalmán E, Gerekliglu C, Yeral M, Buyukkurt N, Aytan P, et al. Quantum cell expansion system: safe and rapid expansion. *Cytotherapy* 2017;19:1246–7.
- [10] Jones M, Varella-Garcia M, Skokan M, Bryce S, Schowinsky J, Peters R, et al. Genetic stability of bone marrow-derived human mesenchymal stromal cells in the Quantum System. *Cytotherapy* 2013;15:1323.
- [11] Mennan C, Garcia J, Roberts S, Hulme C, Wright K. A comprehensive characterisation of large-scale expanded human bone marrow and umbilical cord mesenchymal stem cells. *Stem Cell Res Ther* 2019;10.
- [12] Zhang J, Buller B, Zhang Z, Zhang Y, Lu M, Rosene D, et al. Exosomes derived from bone marrow mesenchymal stromal cells promote remyelination and reduce neuroinflammation in the demyelinating central nervous system. *Exp Neurol* 2021;347:113895.
- [13] Maeda T, Briggs CM, Datar A, Brantner CA, Hanley PJ, Jonas RA, et al. Influence of administration of mesenchymal stromal cell on pediatric oxygenator performance and inflammatory response. *JTCVS Open* 2021;5:99–107.
- [14] Savelli S, Trombi L, D'Alessandro D, Moscato S, Pacini S, Giannotti S, et al. Pooled human serum: a new culture supplement for bioreactor-based cell therapies. Preliminary results. *Cytotherapy* 2018;20:556–63.
- [15] Russell AL, Lefavor RC, Zubair AC. Characterization and cost-benefit analysis of automated bioreactor-expanded mesenchymal stem cells for clinical applications. *Transfusion (Paris)* 2018;58:2374–82.
- [16] Nold P, Brendel C, Neubauer A, Bein G, Hackstein H. Good manufacturing practice-compliant animal-free expansion of human bone marrow derived mesenchymal stroma cells in a closed hollow-fiber-based bioreactor. *Biochem Biophys Res Commun* 2013;430:325–30.
- [17] Lechanteur C, Baila S, Janssens ME, Giet O, Briquet A, Baudoux E, et al. Large-scale clinical expansion of mesenchymal stem cells in the GMP-compliant, closed automated Quantum<sup>®</sup> cell expansion system: comparison with expansion in traditional T-flasks. *J Stem Cell Res Ther* 2014;04.
- [18] Haack-Sørensen M, Juhl M, Follin B, Harary Søndergaard R, Kirchoff M, Kastrup J, et al. Development of large-scale manufacturing of adipose-derived stromal cells for clinical applications using bioreactors and human platelet lysate. *Scand J Clin Lab Invest* 2018;78:293–300.
- [19] Pirrone C, Gobetti A, Caprara C, Bernardini G, Gornati R, Soldati G. Chondrogenic potential of hASCs expanded in flask or in a hollow-fiber bioreactor. *J Stem Cell Res Med* 2017;2.
- [20] Schive SW, Fjukstad R, Josefsen D, Katavić M, Abadpour S, Gullestad H, et al. Automated isolation and expansion of human adipose tissue-derived stem cells for a seamless translation into clinical trials. *CellR4* 2018;6:e2519.
- [21] Coccè V, la Monica S, Bonelli M, Alessandri G, Alfieri R, Lagrasta CA, et al. Inhibition of human malignant pleural mesothelioma growth by mesenchymal stromal cells. *Cells* 2021;10:1427.
- [22] Mizukami A, de Abreu Neto MS, Moreira F, Fernandes-Platzgummer A, Huang YF, Milligan W, et al. A fully-closed and automated hollow fiber bioreactor for clinical-grade manufacturing of human mesenchymal stem/stromal cells. *Stem Cell Rev Rep* 2018;14:141–3.
- [23] Møller-Hansen M, Larsen AC, Toft PB, Lynggaard CD, Schwartz C, Bruunsgaard H, et al. Safety and feasibility of mesenchymal stem cell therapy in patients with aqueous deficient dry eye disease. *Ocular Surface* 2021;19:43–52.
- [24] Bellio MA, Kanashiro-Takeuchi RM, Takeuchi L, Kulandavelu S, Lee Y-S, Balkan W, et al. Systemic delivery of large-scale manufactured Wharton's Jelly mesenchymal stem cell-derived extracellular vesicles improves cardiac function after myocardial infarction. *J Cardiovasc Aging* 2022;2:9.
- [25] Gutova M, Cheng JP, Adhikarla V, Tsururyan L, Barish ME, Rockne RC, et al. Intranasally administered L-Myc-immortalized human neural stem cells migrate to primary and distal sites of damage after cortical impact and enhance spatial learning. *Stem Cells Int* 2021;2021.
- [26] Lambrechts T, Papanthiou I, Rice B, Schrooten J, Luyten FP, Aerts JM. Large-scale progenitor cell expansion for multiple donors in a monitored hollow fibre bioreactor. *Cytotherapy* 2016;18:1219–33.
- [27] Sheu J, Beltzer J, Fury B, Wilczek K, Tobin S, Falconer D, et al. Large-scale production of lentiviral vector in a closed system hollow fiber bioreactor. *Mol Ther Methods Clin Dev* 2015;2:15020.
- [28] Vang B, Frank N, Jones M, Nankervis B, Coeshott C. Expansion and cellular characterization of primary human adherent cells in the Quantum<sup>®</sup> Cell Expansion System, a hollow-fiber bioreactor system. *J Biol Methods* 2020;7:e130.
- [29] Uslu U, Erdmann M, Wiesinger M, Schuler G, Schuler-Thurner B. Automated Good Manufacturing Practice-compliant generation of human monocyte-derived dendritic cells from a complete apheresis product using a hollow-fiber bioreactor system overcomes a major hurdle in the manufacture of dendritic cells for cancer vaccines. *Cytotherapy* 2019;21:1166–78.
- [30] Trivedi A, Miyazawa B, Gibb S, Valanoski K, Vivona L, Lin M, et al. Bone marrow donor selection and characterization of MSCs is critical for pre-clinical and clinical cell dose production. *J Transl Med* 2019;17:1–16.
- [31] Coeshott C, Vang B, Jones M, Nankervis B. Large-scale expansion and characterization of CD3+ T-cells in the Quantum<sup>®</sup> Cell Expansion System. *J Transl Med* 2019;17:1–13.
- [32] Nankervis B, Jones M, Vang B, Brent Rice R, Coeshott C, Beltzer J. Optimizing T cell expansion in a hollow-fiber bioreactor. *Curr Stem Cell Rep* 2018;4:46–51.
- [33] Jones M, Nankervis B, Roballo KS, Pham H, Bushman J, Coeshott C. A comparison of automated perfusion- and manual diffusion-based human regulatory T cell expansion and functionality using a soluble activator complex. *Cell Transplant* 2020;29.
- [34] Schiess M, Suescun J, Doursout MF, Adams C, Green C, Salterelli JG, et al. Allogeneic bone marrow-derived mesenchymal stem cell safety in idiopathic Parkinson's disease. *Movement Disorders* 2021;36:1825–34.
- [35] Mount NM, Ward SJ, Kefalas P, Hyllner J. Cell-based therapy technology classifications and translational challenges. *Philos Trans R Soc B Biol Sci* 2015;370:20150017.
- [36] Pigeau GM, Csaszar E, Dulgar-Tulloch A. Commercial scale manufacturing of allogeneic cell therapy. *Front Med (Lausanne)* 2018;5:233.
- [37] Kabat M, Bobkov I, Kumar S, Grumet M. Trends in mesenchymal stem cell clinical trials 2004-2018: is efficacy optimal in a narrow dose range? *Stem Cells Transl Med* 2020;9:17–27.
- [38] Coopman K, Medcalf N. From production to patient: challenges and approaches for delivering cell therapies. *Stembook* 2008.
- [39] Lechanteur C, Briquet A, Bettonville V, Baudoux E, Beguin Y. MSC manufacturing for academic clinical trials: from a clinical-grade to a full GMP-compliant process. *Cells* 2021;10:1320.
- [40] Mennan C, Garcia J, McCarthy H, Owen S, Perry J, Wright K, et al. Human articular chondrocytes retain their phenotype in sustained hypoxia while normoxia promotes their immunomodulatory potential. *Cartilage* 2019;10:467–79.
- [41] Adesida AB, Mulet-Sierra A, Jomha NM. Hypoxia mediated isolation and expansion enhances the chondrogenic capacity of bone marrow mesenchymal stromal cells. *Stem Cell Res Ther* 2012;3:1–13.
- [42] Meretoja V v, RL Dahlin, Wright S, Kasper FK, Mikos AG. The effect of hypoxia on the chondrogenic differentiation of co-cultured articular chondrocytes and mesenchymal stem cells in scaffolds. *Biomaterials* 2013;34:4266–73.
- [43] Chen W, Zhuo Y, Duan D, Lu M. Effects of hypoxia on differentiation of mesenchymal stem cells. *Curr Stem Cell Res Ther* 2019;15:332–9.
- [44] Rosová I, Dao M, Capoccia B, Link D, Nolta JA. Hypoxic preconditioning results in increased motility and improved therapeutic potential of human mesenchymal stem cells. *Stem Cells* 2008;26:2173–82.
- [45] Parisi L, Toffoli A, Ghezzi B, Mozzoni B, Lumetti S, Macaluso GM. A glance on the role of fibronectin in controlling cell response at biomaterial interface. *Jpn Dent Sci Rev* 2020;56:50–5.
- [46] Schwartz I, Seger D, Shaltiel S. Vitronectin. *Int J Biochem Cell Biol* 1999;31:539–44.
- [47] Prittie J. The role of cryoprecipitate in human and canine transfusion medicine. *J Vet Emerg Crit Care* 2021;31:204–14.
- [48] Harrison P, Hopkins T, Hulme C, McCarthy H, Wright K. Chondrocyte isolation and expansion. *Methods Mol Biol* 2023;2598:9–19.

- [49] Dominici M, le Blanc K, Mueller I, Slaper-Cortenbach I, Marini FC, Krause DS, et al. Minimal criteria for defining multipotent mesenchymal stromal cells. The International Society for Cellular Therapy position statement. *Cytotherapy* 2006; 8:315–7.
- [50] Krampera M, Galipeau J, Shi Y, Tarte K, Sensebe L. Immunological characterization of multipotent mesenchymal stromal cells—The International Society for Cellular Therapy (ISCT) working proposal. *Cytotherapy* 2013;15:1054–61.
- [51] Perry J, Roelofs AJ, Mennan C, McCarthy HS, Richmond A, Clark SM, et al. Human mesenchymal stromal cells enhance cartilage healing in a murine joint surface injury model. *Cells* 2021;10:1999.
- [52] Smith TA. CAR-T cell expansion in a Xuri cell expansion system W25. *Methods Mol Biol* 2020;2086:151–63.
- [53] Wang X, Urak R, Walter M, Guan M, Han T, Vyas V, et al. Large-scale manufacturing and characterization of CMV-CD19CAR T cells. *J Immunother Cancer* 2022;10.
- [54] Mock U, Nickolay L, Philip B, Cheung GWK, Zhan H, Johnston ICD, et al. Automated manufacturing of chimeric antigen receptor T cells for adoptive immunotherapy using CliniMACS prodigy. *Cytotherapy* 2016;18:1002–11.
- [55] Abou-el-Enin M, Elsallab M, Feldman SA, Fesnak AD, Heslop HE, Marks P, et al. Scalable manufacturing of CAR T cells for cancer immunotherapy. *Blood Cancer Discov* 2021;2:408–22.
- [56] Yáñez-Mó M, Siljander PRM, Andreu Z, Zavec AB, Borràs FE, Buzás EI, et al. Biological properties of extracellular vesicles and their physiological functions. *J Extracell Vesicles* 2015;4:1–60.
- [57] Théry C, Amigorena S, Raposo G, Clayton A. Isolation and characterization of exosomes from cell culture supernatants and biological fluids. *Curr Protoc Cell Biol* 2006. Chapter 3: Unit 3.22.
- [58] Williams AM, Dennahy IS, Bhatti UF, Halaweish I, Xiong Y, Chang P, et al. Mesenchymal stem cell-derived exosomes provide neuroprotection and improve long-term neurologic outcomes in a swine model of traumatic brain injury and hemorrhagic shock. *J Neurotrauma* 2019;36:54–60.
- [59] Mendt M, Kamerkar S, Sugimoto H, McAndrews KM, Wu CC, Gagea M, et al. Generation and testing of clinical-grade exosomes for pancreatic cancer. *JCI Insight* 2018;3.
- [60] Go V, Bowley BGE, Pessina MA, Zhang ZG, Chopp M, Finklestein SP, et al. Extracellular vesicles from mesenchymal stem cells reduce microglial-mediated neuroinflammation after cortical injury in aged Rhesus monkeys. *Geroscience* 2020; 42:1–17.
- [61] Moore TL, Bowley BGE, Pessina MA, Calderazzo SM, Medalla M, Go V, et al. Mesenchymal derived exosomes enhance recovery of motor function in a monkey model of cortical injury. *Restor Neurol Neurosci* 2019;37:347–62.
- [62] Gudbergsson JM, Johnsen KB, Skov MN, Duroux M. Systematic review of factors influencing extracellular vesicle yield from cell cultures. *Cytotechnology* 2016;68:579–92.
- [63] Lötvall J, Hill AF, Hochberg F, Buzás EI, Vizio D di, Gardiner C, et al. Minimal experimental requirements for definition of extracellular vesicles and their functions: a position statement from the International Society for Extracellular Vesicles. *J Extracell Vesicles* 2014;3.
- [64] Théry C, Witwer KW, Aikawa E, Alcaraz MJ, Anderson JD, Andriantsitohaina R, et al. Minimal information for studies of extracellular vesicles 2018 (MISEV2018): a position statement of the International Society for Extracellular Vesicles and update of the MISEV2014 guidelines. *J Extracell Vesicles* 2018;7.
- [65] Portnow J, Synold TW, Badie B, Tirughana R, Lacey SF, D'Apuzzo M, et al. Neural stem cell-based anticancer gene therapy: a first-in-human study in recurrent high-grade glioma patients. *Clin Cancer Res* 2017;23:2951–60.
- [66] Portnow J, Badie B, Suzette Blanchard M, Kilpatrick J, Tirughana R, Metz M, et al. Feasibility of intracerebrally administering multiple doses of genetically modified neural stem cells to locally produce chemotherapy in glioma patients. *Cancer Gene Ther* 2021;28:294.
- [67] Kastrup J, Haack-Sørensen M, Juhl M, Harary Søndergaard R, Follin B, Drozd Lund L, et al. Cryopreserved off-the-shelf allogeneic adipose-derived stromal cells for therapy in patients with ischemic heart disease and heart failure—a safety study. *Stem Cells Transl Med* 2017;6:1963–71.
- [68] Kastrup J, Schou M, Gustafsson I, Nielsen OW, Møgelvang R, Kofoed KF, et al. Rationale and design of the first double-blind, placebo-controlled trial with allogeneic adipose tissue-derived stromal cell therapy in patients with ischemic heart failure: a phase II Danish multicentre study. *Stem Cells Int* 2017;2017.
- [69] Paitazoglou C, Bergmann MW, Vrtovec B, Chamuleau SAJ, van Klarenbosch B, Wojakowski W, et al. Rationale and design of the European multicentre study on Stem Cell therapy in Ischemic Non-treatable Cardiac diseases (SCIENCE). *Eur J Heart Fail* 2019;21:1032.
- [70] Hess DC, Sila CA, Furlan AJ, Wechsler LR, Switzer JA, Mays RW. A double-blind placebo-controlled clinical evaluation of MultiStem for the treatment of ischemic stroke. *Int J Stroke* 2013;9:381–6.
- [71] Hess DC, Wechsler LR, Clark WM, Savitz SI, Ford GA, Chiu D, et al. Safety and efficacy of multipotent adult progenitor cells in acute ischaemic stroke (MASTERS): a randomised, double-blind, placebo-controlled, phase 2 trial. *Lancet Neurol* 2017;16:360–8.
- [72] Osanai T, Houkin K, Uchiyama S, Minematsu K, Taguchi A, Terasaka S. Treatment evaluation of acute stroke for using in regenerative cell elements (TREASURE) trial: Rationale and design. *Int J Stroke* 2018;13:444–8.
- [73] Kincaid R. Top Line Results of the ONE-BRIDGE Study in Patients with ARDS. Tokyo: 2021. <https://moneyworld.jp/disci-pdf/tdnet/2021080548046701GENERAL.pdf> [accessed 01.08.22].
- [74] Bellingan G, Jacono F, Bannard-Smith J, Brealey D, Meyer N, Thickett D, et al. Safety and efficacy of multipotent adult progenitor cells in acute respiratory distress syndrome (MUST-ARDS): a multicentre, randomised, double-blind, placebo-controlled phase 1/2 trial. *Intensive Care Med* 2022;48:36–44.
- [75] Marino D, Ronfard V. The future of bioengineering personalised skin on demand for large skin defects. *J Wound Care* 2021;30:878–9.